



Cole-Parmer®

MIXING SOLUTIONS

IKA® RW 20 Digital Dual-Range Mixer

~~List Price : Rs. 31,600/-*~~

Special Price : Rs. 16,900/-*

1 Mixer

MK-50705-10

- Through-shaft design provides virtual limitless adjustment for impeller depth
- Two speed ranges to suit a variety of mixing applications
- Maintenance-free motor rated for continuous duty

Use this versatile, affordable mixer for liquid volumes up to 20 litres and viscosities up to 10,000 cps (approximately the viscosity of honey). Select speed from 60 to 2000 rpm in two ranges for a variety of applications. Through-shaft design permits adjustment of shaft length and impeller depth—allows you to place or remove mixing vessels without repositioning the mixer.

The motor and drive unit are conveniently integrated into a compact aluminum and thermoplastic housing with easy-to-use controls. The quiet, fan-cooled motor features thermal overload shut-off and an efficient heat dissipating design to ensure a long service life. Motor is maintenance-free and rated for continuous duty. Mixer features stepless speed control. The 0.5" LCD shows rotational speed in rpm and enables you to precisely control and repeat process parameters.

What's included: chuck (accepts up to 0.375" [10-mm] dia shafts), 6.25"L x 0.5" dia mounting rod and a 6.5-ft, three wire power cord.



Specifications

Speed Range 1	: 60 to 500 rpm
Speed Range 2	: 240 to 2000 rpm
Motor Power Input	: 70 Watts
Dimensions	: 3.5"W x 11.5"H x 7.375"D
HP	: 1/10
Max torque	: 259 in-oz
Power	: 230 VAC

Ready to use Package 1 + 2 + 3 + 4

2 Four-Blade Propeller

Standard stirring element for drawing the material to be mixed from the top to the bottom. Local shearing forces. Generates axial flow in the vessel. Use at medium to high speeds.

Specifications

Blade diameter	: 2" (50 mm)
Shaft diameter	: 0.3125" (8 mm)
Shaft length	: 13.75" (350 mm)
Speed max	: 2000 rpm 1/min



~~List Price : Rs. 60,800/-*~~

Special Price : Rs. 22,000/-*

MK-15940-64

CE

4 Boss Head Clamp

Specifications

Clamp range	: 0.25" to 0.625" (6 to 16 mm)
Material	: Cast aluminum




3 Plate Stand

With slip-resistant foil.

Specifications

Support rod dia	: 0.625" (16 mm)
Max load	: 11 lb (5 kg)
Dimensions	: 0.75" x 1.25" x 32"
(W X D x H)	(200 x 316 x 800 mm)

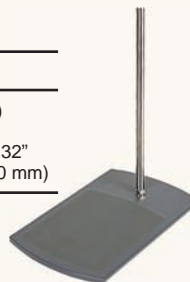

*Terms & Conditions: Prices quoted are Ex-works warehouse • Taxes & other levies extra as applicable • Offer valid only upto 31st March 2010

Cole-Parmer India (Div. of Scientific Products India Pvt Ltd)

Godrej Coliseum 101 A - 101 B, Somaiya Hospital Road, Sion (E), Mumbai - 400022, India

Tel: +91-22-67162222 • Fax: +91-22-67162211 • Email: info@coleparmer.in • Web: www.ColeParmer.in